

Application No. 09/845,985

REMARKS

Claims 1-3, 6-9, 12, 14-21, 48-50, 52-56 and 58-61 are pending. Claims 1 and 21 are amended to correct an obvious typographical error. All of the pending claims stand rejected. Applicant respectfully requests reconsideration of the rejection based on the following comments.

Rejection Over Kamauchi and Manev

The Examiner rejected claims 1-3, 6, 7, 12, 14-17, 19-21, 48-50, 52, 53, 55, 56 and 58-61 as being unpatentable under 35 U.S.C. § 103(a) over U.S. patent 5,538,814 to Kamauchi (Kamauchi) in view of U.S. patent 5,789,115 to Manev (Manev). Except for insignificant changes to the notation in references to the cited art, Applicant notes that the text of this rejection seems identical to the rejection filed on December 13, 2006. Applicant incorporates by reference their arguments presented in Preliminary Remarks dated September 14, 2009 as well as their Appeal Brief dated May 14, 2007 with two corrected pages filed on June 26, 2007, Reply Brief filed November 29, 2007 and Request for Rehearing filed on December 29, 2009. Here, Applicant focuses on the Examiner's response on pages 2 and 3 of the Office Action relating to the new Horne Declaration.

The Examiner has raised questions regarding the analysis presented in the Declaration without any documentary support whatsoever. With all due respect, the Examiner's position is not well founded. The Examiner calls into question to issues relating to the mass of the particles and asserts that the claimed invention is not commensurate with the evidence presented because mass is not in the claim. The mass is only involved in the data as related to relating the calculated volume distributions with number distributions. The volume is scaled based on the density to obtain the mass. As explained in the Declaration, the mass percent in the Horne Declaration is used as a substitute for volume. The Office Action asserts that "the mass does not

Application No. 09/845,985

account for any discrepancies in shapes or the porosity of the material but assumes a perfect material of the same density throughout and a perfect sphere."

With respect to the porosity, this is a known property. The micrographs are examined to determine porosity. No porosity is visible in the micrographs. Therefore, porosity is ruled out as a significant issue. With respect to density, the known composition and lack of porosity imply that density of the material is consistent for the particles.

With respect to the non-spherical nature of the particles, this is an issue, but this issue supports Applicant's position significantly. Since light scattering measures a hydration sphere, the technique is directed to spheres. Applicant has attached a Technical Note from Malvern Instruments on Dynamic Light Scattering. The issue of non-spherical particles is explicitly discussed. Since an examination of the micrographs indicates that the particles for the most part are roughly spherical, an even more significant issue is the measurement of a hydrodynamic diameter rather than a physical diameter.

As discussed in Applicant's specification, the particle size, as specified in paragraph [0119], the particle sizes are determined by transmission electron microscopy, which is the most accurate technique. Also, for non-spherical particles, the average is taken along the principle axes of the particle.

The Horne Declaration of September 22, 2009 directly addresses this issue in section 7. There it is stated that "In my experience, light scattering will result in a narrower distribution than a TEM analysis for at least two reasons. First, there is small particle drop out in the light scattering measurements, so the small particle part of the distribution is under represented or absent, [sic] Also, agglomeration tends to involve more smaller particles so that you bias the distribution toward the middle." Thus, light scattering measurements, with its imprecision relating to particle shapes, tends to give a significantly narrow particle size distribution than the proper transmission electron micrograph (TEM) approach. The particle shape issue will also

Application No. 09/845,985

tend to produce a narrow distribution due to bias toward larger particle sizes in the measurements. But even with these biases toward a narrower distribution, the dynamic light scattering measurements still produced a particle size distribution that was broader than the claimed distributions.

Since the issues of imprecision of dynamic light scattering generally were addressed in the Horne Declaration, the conclusions of the Horne Declaration are clearly sound. The Examiner has failed to raise any significant issues to rebut the conclusions of the Horne Declaration. Therefore, it is clear that the teachings of the cited references fail to establish a reasonable expectation of success with respect to Applicant's claimed invention.

The Office Action states that regarding "the particle sizes, the declaration illustrates that the sizes taught by Kamauchi are obtainable and therefore a reasonable expectation of success does exist. As taught by the prior arts of Kamauchi and Manev, longer grinding times will produce smaller particles and thus lower the average particle size. Applicant has not provided evidence or reasoning why this teaching would not produce the claimed sizes. Furthermore, applicant has not provided evidence why Kamauchi's example 13 with an average particle size of 0.01 microns would not meet the claimed invention."

With all due respect, this analysis completely ignores features of the claimed invention, which have been discussed in great detail over extended periods of prosecution. All of the pending claims have features relating to the uniformity of the particles. The uniformity feature is a very important characteristic of the claimed materials. This is discussed in detail in the briefs and the earlier file history. For example, claims 1 and 21 indicate that essentially no particles have a diameter greater than about 5 times the average particle size. Also, claims 55 and 58 indicate uniformity with respect to a distribution of particle sizes such that at least about 95 percent of the particles have a diameter greater than about 40 percent of the average diameter and less than about 160 percent of the average diameter. The uniformity of particle size distribution is an

Application No. 09/845,985

independent but related parameter of the materials. For example, a soccer ball and a ping pong ball may have the same average diameter as two softballs, but you cannot play softball with the soccer ball and ping pong ball. There are different articles, and the powders with different particle size distributions are different composition of matter with different properties. With all due respect, Applicant has presented considerable un-refuted evidence that the teachings of the art do not place the claimed materials in the public domain. Therefore, it is clear that the combined teachings of the cited references do not render Applicant's claimed invention obvious.

In view of the comments above, Applicant respectfully requests withdrawal of the rejection of 1-3, 6, 7, 12, 14-17, 19-21, 48-50, 52, 53, 55, 56 and 58-61 as being unpatentable under 35 U.S.C. § 103(a) over U.S. patent 5,538,814 to Kamauchi (Kamauchi) in view of U.S. patent 5,789,115 to Manev (Manev).

Rejection Over Goodenough, Kamauchi and Manev

The Examiner rejected claims 8, 9 and 18 under 35 U.S.C. § 103(a) as being unpatentable over U.S. patent 5,910,382 to Goodenough (Goodenough) in view of Kamauchi and further in view of Manev. Except for insignificant changes to the notation of references to the cited art, Applicant notes that the text of this rejection seems identical to the rejection filed on December 13, 2006. Applicant incorporates by reference their arguments presented in Preliminary Remarks dated September 14, 2009 as well as their Appeal Brief dated May 14, 2007 with two corrected pages filed on June 26, 2007, Reply Brief filed November 29, 2007 and Request for Rehearing filed on December 29, 2009. Here, Applicant focuses on the Examiner's response on pages 2 and 3 of the Office Action relating to the new Horne Declaration. These issues regarding the Horne Declaration have been addressed in detail above. In view of the comments above and the analysis incorporated herein by reference, the combined teachings of the cited references clearly to not render Applicant's claimed invention obvious. Applicant respectfully requests withdrawal

Application No. 09/845,985

of the rejection of claims 8, 9 and 18 under 35 U.S.C. § 103(a) as being unpatentable over Goodenough in view of Kamauchi and further in view of Manev.

Rejection Of Claim 54

The Examiner rejected claim 54 under 35 U.S.C. § 103(a) as being unpatentable over Kamauchi in view of Manev and further in view of U.S. patent 5,232,794 to Krumpelt et al. (Krumpelt). The Examiner cited Krumpelt for teaching the use of AlPO_4 "for the conduction of lithium in batteries." Claim 54 depends from claim 21. First, Krumpelt is silent with respect to the formation of highly uniform nanoparticles or appropriate ways of forming such highly uniform nanoparticles. Thus, Krumpelt does not make up for the deficiencies of the combined teachings of Kamauchi and Manev.

Furthermore, ion conduction, as taught in Krumpelt, relates to the formation of an electrolyte, see column 7, lines 20-30. Solid electrolytes in batteries are replacements for separators. See, for example, US 2009/0029264 to Nakazawa et al, entitled "Thin Film Solid Secondary Cell. Separators must provide electrical insulation and a physical barrier between the anode and cathode. As described at the top of column 6 of Krumpelt, the pellets are pressed and sintered. Thus, Krumpelt teaches away from nanoparticles. (KSR)

The combined teachings of Kamauchi, Manev and Krumpelt clear do not render claim 54 *prima facie* obvious. Applicant respectfully requests withdrawal of the rejection of claim 54 under 35 U.S.C. § 103(a) as being unpatentable over Kamauchi in view of Manev and further in view of Krumpelt.


CONCLUSIONS

In view of the foregoing, it is submitted that this application is in condition for allowance. Favorable consideration and prompt allowance of the application are respectfully requested.

Application No. 09/845,985

The Examiner is invited to telephone the undersigned if the Examiner believes it would be useful to advance prosecution.

Respectfully submitted,

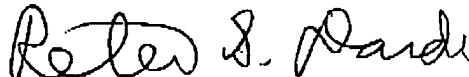A handwritten signature in black ink, appearing to read "Peter S. Dardi". The signature is fluid and cursive, with the first name "Peter" and last name "Dardi" clearly distinguishable.

Peter S. Dardi, Ph.D.
Registration No. 39,650

Customer No. 62274
Dardi & Herbert, PLLC
US Bank Plaza, Suite 2000
220 South 6th Street
Minneapolis, Minnesota 55402
Telephone: (404) 949-5730



technical

note

Dynamic Light Scattering: An Introduction in 30 Minutes



Introduction

Dynamic Light Scattering (sometimes referred to as Photon Correlation Spectroscopy or Quasi-Elastic Light Scattering) is a technique for measuring the size of particles typically in the sub micron region.

Brownian Motion

DLS measures Brownian motion and relates this to the size of the particles. Brownian motion is the random movement of particles due to the bombardment by the solvent molecules that surround them. Normally DLS is concerned with measurement of particles suspended within a liquid.

The larger the particle, the slower the Brownian motion will be. Smaller particles are "kicked" further by the solvent molecules and move more rapidly. An accurately known temperature is necessary for DLS because knowledge of the viscosity is required (because the viscosity of a liquid is related to its temperature). The temperature also needs to be stable, otherwise convection currents in the sample will cause non-random movements that will ruin the correct interpretation of size.

The velocity of the Brownian motion is defined by a property known as the translational diffusion coefficient (usually given the symbol, D).

The Hydrodynamic Diameter

The size of a particle is calculated from the translational diffusion

coefficient by using the Stokes-Einstein equation:

$$d(H) = \frac{kT}{3\pi\eta D}$$

where:-

$d(H)$ = hydrodynamic diameter

D = translational diffusion coefficient

k = Boltzmann's constant

T = absolute temperature

η = viscosity

Note that the diameter that is measured in DLS is a value that refers to how a particle diffuses within a fluid so it is referred to as a hydrodynamic diameter. The diameter that is obtained by this technique is the diameter of a sphere that has the same translational diffusion coefficient as the particle.

The translational diffusion coefficient will depend not only on the size of the particle "core", but also on any surface structure, as well as the concentration and type of ions in the medium. Factors that affect the diffusion speed of particles are discussed in the following sections.

Ionic Strength of Medium

The ions in the medium and the total ionic concentration can affect the particle diffusion speed by changing the thickness of the electric double layer called the Debye length (K^{-1}). Thus a low conductivity medium will produce an extended double layer of ions around the particle,

reducing the diffusion speed and resulting in a larger, apparent hydrodynamic diameter. Conversely, higher conductivity media will suppress the electrical double layer and the measured hydrodynamic diameter.

The performance of a DLS instrument is normally verified by measurement of a suitable polystyrene latex standard. If the standard needs to be diluted prior to measurement, then dilution in an appropriate medium is important. The International Standard on DLS (ISO13321 Part 8 1996) says that dilution of any polystyrene standard should be made in 10mM NaCl. This concentration of salt will suppress the electrical double layer and ensure that the hydrodynamic diameter reported will be the same as the hydrodynamic diameter on the certificate or the expected diameter.

Surface Structure

Any change to the surface of a particle that affects the diffusion speed will correspondingly change the apparent size of the particle. An adsorbed polymer layer projecting out into the medium will reduce the diffusion speed more than if the polymer is lying flat on the surface. The nature of the surface and the polymer, as well as the ionic concentration of the medium can affect the polymer conformation, which in turn can change the apparent size by several nanometres.

Non-Spherical Particles

All particle-sizing techniques have an inherent problem in describing the size of non-spherical particles. The sphere is the only object whose size





technical note

can be unambiguously described by a single figure.

Different techniques are sensitive to different properties of the particle, e.g. projected area, density, scattering intensity, and in general will produce different mean sizes and size distributions for any given sample. Even the size in a microscope image will depend on parameters set such as edge contrast etc. It is important to understand that none of these results are inherently "correct".

The hydrodynamic diameter of a non-spherical particle is the diameter of a sphere that has the same translational diffusion speed as the particle.

If the shape of a particle changes in a way that affects the diffusion speed, then the hydrodynamic size will change. For example, small changes in the length of a rod-shaped particle will directly affect the size, whereas changes in the rod's diameter, which will hardly affect the diffusion speed, will be difficult to detect.

The conformation of proteins and macromolecules are usually dependent on the exact nature of the dispersing medium. As conformational changes will usually affect the diffusion speed, DLS is a very sensitive technique for detecting these changes.

Light Scattering Theories

Rayleigh Scattering

If the particles are small compared to the wavelength of the laser used (typically less than $d = \lambda/10$ or around 60nm for a He-Ne laser), then the scattering from a particle illuminated by a vertically polarised laser will be essentially isotropic, i.e. equal in all directions.

The Rayleigh approximation tells us that $I \propto d^6$ and also that $I \propto 1/\lambda^4$, where I = intensity of light scattered, d = particle diameter and λ = laser

wavelength. The d^6 term tells us that a 50nm particle will scatter 10^6 or one million times as much light as a 5nm particle. Hence there is a danger that the light from the larger particles will swamp the scattered light from the smaller ones. This d^6 factor also means it is difficult with DLS to measure, say, a mixture of 1000nm and 10nm particles because the contribution to the total light scattered by the small particles will be extremely small. The inverse relationship to λ^4 means that a higher scattering intensity is obtained as the wavelength of the laser used decreases.

Mie Theory

When the size of the particles becomes roughly equivalent to the wavelength of the illuminating light, then a complex function of maxima and minima with respect to angle is

observed.

Figure 1 shows the theoretical plot of the log of the relative scattering intensity versus particle size at angles of 173° (the detection angle of the Zetasizer Nano S and Nano ZS in aqueous media) and 90° (the detection angle of the Nano S90 and Nano ZS90) assuming a laser wavelength of 633nm, real refractive index of 1.59 and an imaginary refractive index of 0.001. Mie theory is the only theory that explains correctly the maxima and minima in the plot of intensity with angle and will give the correct answer over all wavelengths, sizes and angles. Mie theory is used in the Nano software for conversion of the intensity distribution into volume.

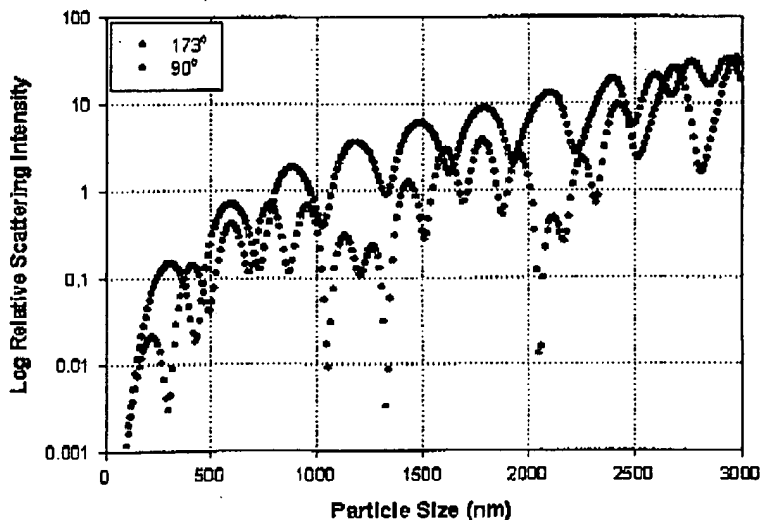


Figure 1: Theoretical plot of the log of the relative intensity of scattering versus particle size at angles of 173° (the detection angle of the Nano S, and Nano ZS in aqueous media) and 90° (the detection angle of the Nano S90 and Nano ZS90) assuming a laser beam at a wavelength of 633nm, real refractive index of 1.59 and an imaginary refractive index of 0.001



technical note

How DLS Works

In dynamic light scattering, the speed at which the particles are diffusing due to Brownian motion is measured. This is done by measuring the rate at which the intensity of the scattered light fluctuates when detected using a suitable optical arrangement. How do these fluctuations in the intensity of scattered light arise?

Imagine if a cuvette, containing particles which are stationary, is illuminated by a laser and a frosted glass screen is used to view the sample cell. A classical speckle pattern would be seen (figure 2). The speckle pattern will be stationary both in speckle size and position because the whole system is stationary. The dark spaces are where the phase additions of the scattered light are mutually destructive and cancel each other out (figure 3A). The bright blobs of light in the speckle pattern are where the light scattered from the particles arrives with the same phase and interfere constructively to form a bright patch (figure 3B).

For a system of particles undergoing Brownian motion, a speckle pattern is observed where the position of each speckle is seen to be in constant motion. This is because the phase addition from the moving particles is constantly evolving and forming new patterns. The rate at which these intensity fluctuations occur will depend on the size of the particles. Figure 4 schematically illustrates typical intensity fluctuations arising from a dispersion of large particles and a dispersion of small particles. The small particles cause the intensity to fluctuate more rapidly than the large ones.

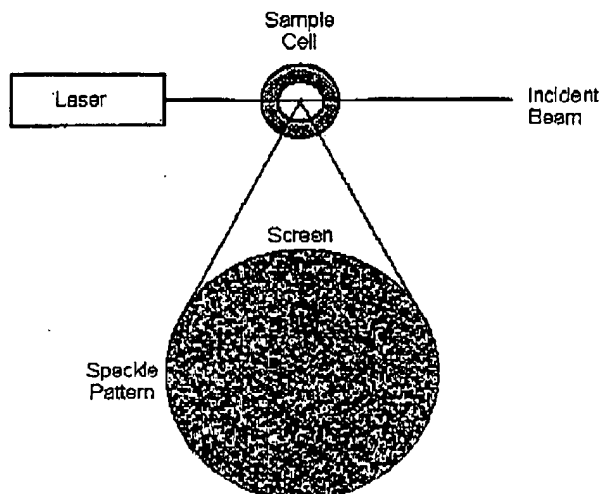


Figure 2: Schematic representation of a speckle pattern

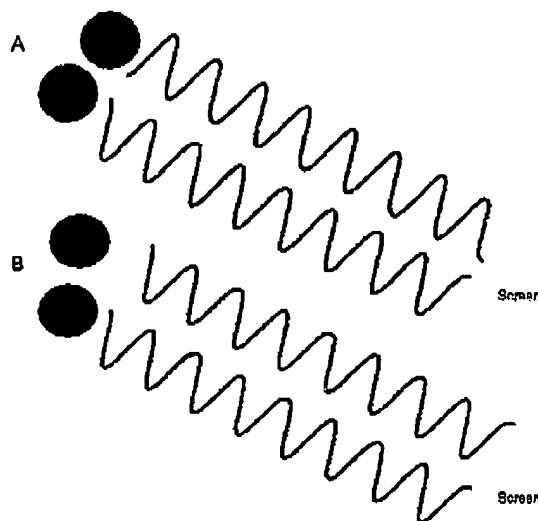


Figure 3: The observed signal depends on the phase addition of the scattered light falling on the detector. In example A, two beams interfere and 'cancel each other out' resulting in a decreased intensity detected. In example B, two beams interfere and 'enhance each other' resulting in an increased intensity detected.



technical note

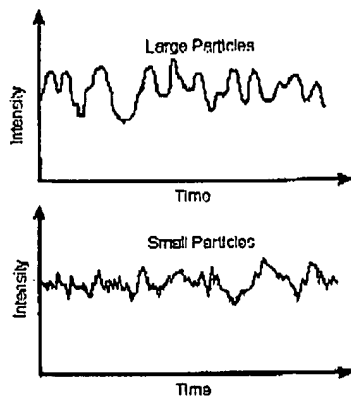


Figure 4: Typical intensity fluctuations for large and small particles

It is possible to directly measure the spectrum of frequencies contained in the intensity fluctuations arising from the Brownian motion of particles, but it is inefficient to do so. The best way is to use a device called a digital auto correlator.

How a Correlator Works

A correlator is basically a signal comparator. It is designed to measure the degree of similarity between two signals, or one signal with itself at varying time intervals.

If the intensity of a signal is compared with itself at a particular point in time and a time much later, then for a randomly fluctuating signal it is obvious that the intensities are not going to be related in any way, i.e. there will be no correlation between the two signals (figure 5). Knowledge of the initial signal intensity will not allow the signal intensity at time $t = \infty$ to be predicted. This will be true of any random process such as diffusion.

However, if the intensity of signal at time t is compared to the intensity a very small time later ($t + \delta t$), there will be a strong relationship or correlation

between the intensities of two signals. The two signals are strongly or well correlated.

If the signal, derived from a random process such as Brownian motion, at t is compared to the signal at $t + 2\delta t$, there will still be a reasonable comparison or correlation between the two signals, but it will not be as good as the comparison at t and $t + \delta t$. The correlation is reducing with time. The period of time δt is usually very small, maybe nanoseconds or microseconds and is called the sample time of the correlator. $t = \infty$ maybe of the order of a millisecond or tens of milliseconds.

If the signal intensity at t is compared with itself then there is perfect correlation as the signals are identical. Perfect correlation is indicated by unity (1.00) and no correlation is indicated by zero (0.00).

If the signals at $t + 2\delta t$, $t + 3\delta t$, $t + 4\delta t$ etc. are compared with the signal at t , the correlation of a signal arriving from a random source will decrease with time until at some time, effectively $t = \infty$, there will be no correlation.

If the particles are large the signal will be changing slowly and the correlation will persist for a long time (figure 6). If the particles are small and moving rapidly then correlation will reduce more quickly (figure 7).


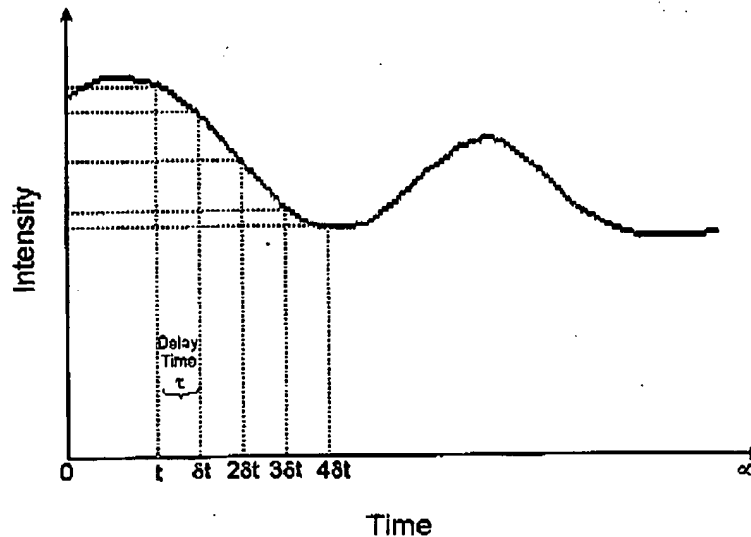


Figure 5: Schematic showing the fluctuation in the intensity of scattered light as a function of time

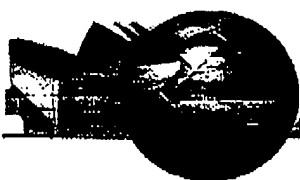

technical

note

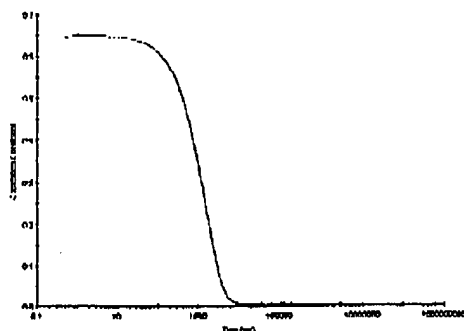


Figure 6: Typical correlogram from a sample containing large particles in which the correlation of the signal takes a long time to decay

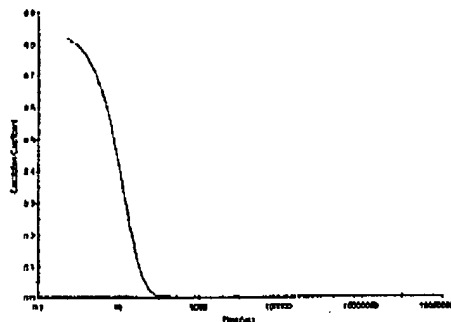


Figure 7: Typical correlogram from a sample containing small particles in which the correlation of the signal decays more rapidly

Viewing the correlogram from a measurement can give a lot of information about the sample. The time at which the correlation starts to significantly decay is an indication of the mean size of the sample. The steeper the line, the more monodisperse the sample is. Conversely, the more extended the decay becomes, the greater the sample polydispersity.

The Correlation Function

It has been seen that particles in a dispersion are in a constant, random Brownian motion and that this causes the intensity of scattered light to fluctuate as a function of time. The correlator used in a PCS Instrument will construct the correlation function $G(\tau)$ of the scattered intensity:

$$G(\tau) = \langle I(t) \cdot I(t+\tau) \rangle$$

Where τ = the time difference (the sample time) of the correlator.

For a large number of monodisperse particles in Brownian motion, the correlation function (given the symbol $G(\tau)$) is an exponential decaying

function of the correlator time delay τ :

$$G(\tau) = A[1 + B \exp(-2\Gamma\tau)]$$

where A = the baseline of the correlation function, B = intercept of the correlation function.

$$\Gamma = Dq^2$$

where D = translational diffusion coefficient

$$q = (4\pi n / \lambda_0) \sin(\theta/2)$$

where n = refractive index of dispersant, λ_0 = wavelength of the laser, θ = scattering angle.

For polydisperse samples, the equation can be written as:

$$G(\tau) = A[1 + B g_1(\tau)^2]$$

where $g_1(\tau)$ = is the sum of all the exponential decays contained in the correlation function.

Obtaining Size Information From the Correlation Function

Size is obtained from the correlation function by using various algorithms. There are two approaches that can be taken (1) fit a single exponential to the correlation function to obtain the

mean size (z-average diameter) and an estimate of the width of the distribution (polydispersity index) (this is called the Cumulants analysis and is defined in ISO13321 Part 8), or (2) fit a multiple exponential to the correlation function to obtain the distribution of particle sizes (such as Non-negative least squares (NNLS) or CONTIN.

The size distribution obtained is a plot of the relative intensity of light scattered by particles in various size classes and is therefore known as an intensity size distribution.

If the distribution by intensity is a single fairly smooth peak, then there is little point in doing the conversion to a volume distribution using the Mie theory. If the optical parameters are correct, this will just provide a slightly different shaped peak. However, if the plot shows a substantial tail, or more than one peak, then Mie theory can make use of the input parameter of sample refractive index to convert the intensity distribution to a volume distribution. This will then give a more realistic view of the importance of the tail or second peak present. In general terms it will be seen that:-

$$d(\text{intensity}) > d(\text{volume}) > d(\text{number})$$





technical note

A very simple way of describing the difference between intensity, volume and number distributions is to consider 2 populations of spherical particles of diameter 5nm and 50nm present in equal numbers (figure 8). If a number distribution of these 2 particle populations is plotted, a plot consisting of 2 peaks (positioned at 5 and 50nm) of a 1 to 1 ratio would be obtained. If this number distribution was converted into volume, then the 2 peaks would change to a 1:1000 ratio (because the volume of a sphere is equal to $4/3\pi(d/2)^3$). If this was further converted into an intensity distribution, a 1:1000000 ratio between the 2 peaks would be obtained (because the intensity of scattering is proportional to d^6 (from Rayleighs approximation)). Remember that in DLS, the distribution obtained from a measurement is based on intensity.

Optical Configuration of a Dynamic Light Scattering Instrument

A typical dynamic light scattering system comprises of six main components. Firstly, a laser \odot provides a light source to illuminate the sample contained in a cell \ominus . For dilute concentrations, most of the laser beam passes through the sample, but some is scattered by the particles within the sample at all angles. A detector \oplus is used to measure the scattered light. In the Zetasizer Nano series, the detector position will be at either 173° or 90° , depending upon the particular model.

The intensity of scattered light must be within a specific range for the detector to successfully measure it. If too much light is detected, then the detector will become saturated. To overcome this, an attenuator \ominus is used to reduce the intensity of the laser source and hence reduce the intensity of scattering. For samples that do not scatter much light,

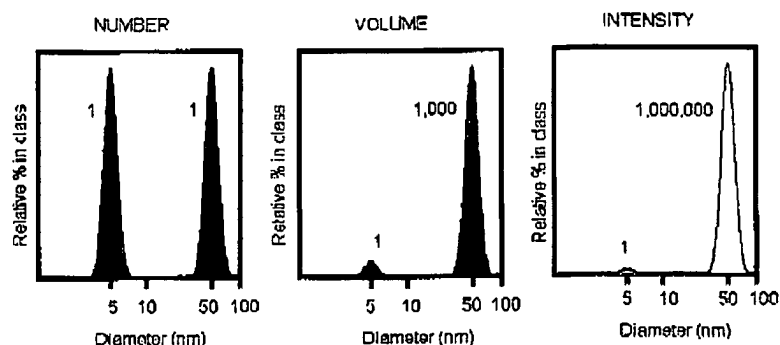


Figure 8: Number, volume and intensity distributions of a bimodal mixture of 5 and 50nm lattices present in equal numbers

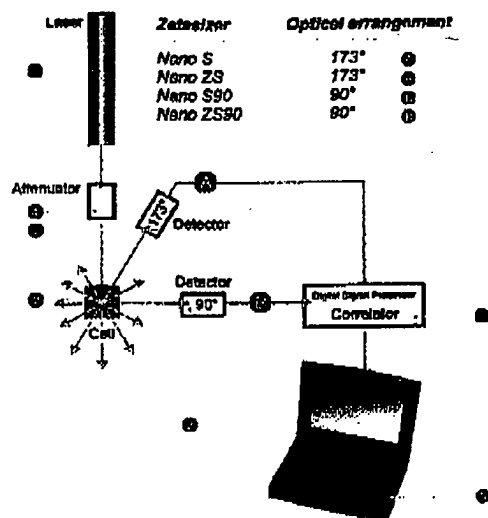


Figure 9: Optical configurations of the Zetasizer Nano series for dynamic light scattering measurements

such as very small particles or samples of low concentration, the amount of scattered light must be increased. In this situation, the attenuator will allow more laser light through to the sample.

For samples that scatter more light, such as large particles or samples at higher concentration, the intensity of scattered light must be decreased. The appropriate attenuator position is automatically determined by the Nano software and covers a transmission range of 100% to 0.0003%.



technical

note

The scattering intensity signal from the detector is passed to a digital processing board called a correlator ③. The correlator compares the scattering intensity at successive time intervals to derive the rate at which the intensity is varying. This correlator information is then passed to a computer ④, where the Nano software will analyze the data and derive size information.

Unique Features of the Zetasizer Nano

Non-Invasive Backscatter Detection (NIBS)

The Nano S and Nano ZS instruments detect the scattering information at 173° . This is known as backscatter detection. In addition, the optics are not in contact with the sample and hence the detection optics are said to be non-invasive. There are several advantages in using non-invasive backscatter detection:

- The laser does not have to travel through the entire sample. This reduces an effect called multiple scattering, where light from one particle is itself scattered by other particles. As the light passes through a shorter path length of the sample, then higher concentrations of sample can be measured.
- Contaminants such as dust particles within the dispersant are typically large compared to the sample size. Large particles mainly scatter in the forward direction. Therefore, by using backscatter detection, the effects of dust are greatly reduced.

Variable Measurement Position For Sizing

The measurement position within the cuvette of the Nano S and Nano ZS can be changed. This measurement position is changed by moving the focusing lens and is determined automatically by the Nano software (figure 10).

For small particles, or samples at low concentrations, it is beneficial to maximise the amount of scattering from the sample. As the laser passes through the wall of the cuvette and into the dispersant, the laser will cause "flare".

This flare may swamp the signal from the scattering particles. Moving the

measurement point away from the cuvette wall towards the centre of the cuvette will remove this effect (figure 10a).

Large particles or samples at high concentrations scatter much more light. In this situation, measuring closer to the cuvette wall will reduce the effect of multiple scattering by minimising the path length over which the scattered light has to pass (figure 10b). The measurement position is determined automatically through a combination of the Intercept of the correlation function and the intensity of the light scattered.

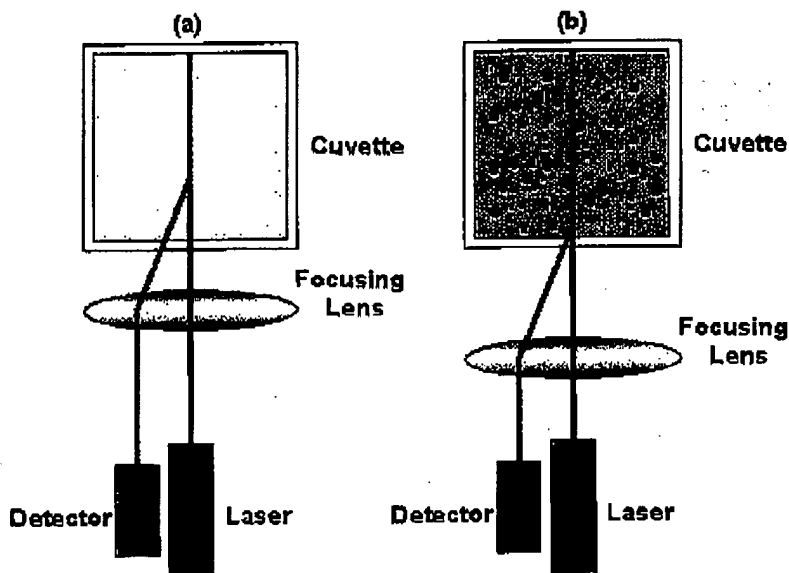


Figure 10: Schematic diagram showing the measurement position for (a) small, weakly scattering samples and for (b) concentrated, opaque samples. The change in measurement position is achieved by moving the focusing lens accordingly

**technical**

note

Additional Reading

[1] International Standard ISO13321
Methods for Determination of Particle
Size Distribution Part 8: Photon
Correlation Spectroscopy,
International Organisation for
Standardisation (ISO) 1998.

[2] Dahneke, B.E. (ed) Measurement
of Suspended Particles by Quasi-
elastic Light Scattering, Wiley, 1983.

[3] Pecora, R. Dynamic Light
Scattering: Applications of Photon
Correlation Spectroscopy, Plenum
Press, 1985.

[4] Washington, C. Particle Size
Analysis In Pharmaceuticals And Other
Industries: Theory And Practice, Ellis
Horwood, England, 1992.

[5] Johnson, C.S. Jr. and Gabriel,
D.A. Laser Light Scattering, Dover
Publications, Inc., New York 1981

Malvern Instruments Ltd

Enigma Business Park • Grovewood Road • Malvern • Worcestershire • UK • WR14 1XZ
Tel: +44 (0)1684 892456 • Fax: +44 (0)1684 892789

Malvern Instruments Worldwide

Sales and service centers in over 50 countries for details visit www.malvern.co.uk/contact

more information at www.malvern.co.uk

8

DLS technical note MK658-01

